



Duckworth & Kent

Ophthalmic Titanium Surgical Instruments


Wallace LRI Set

(Wallace Limbal Relaxing Kit for Reduction of Astigmatism)

Wallace Limbal Relaxing LRI "Kit" allows separate sterilizing capability in order that instruments may be used at the time when limbal relaxing incisions are clinically indicated.

4-620

Wallace LRI Diamond Knife

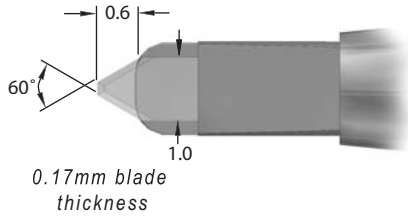

- 0.2mm front flat at tip
- 600 micron setting / extension
- Retractable long handle
- Length 120.0mm

Trifacet blade with single footplate allows easier visibility of knife as it passes through corneal tissue; 600 micron setting. Thin streamlined handle for finger twirling as blade follows arcuate pattern of the limbus.

9-705R

Wallace-Mendez Degree Gauge




- Measures 0° - 180° in 10° increments x 2
- Round handle, length 105.0mm

10° interval marks. Designed for surgeons familiar with use of a standard phoropter.

2-138D

DK Colibri Forceps

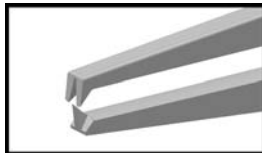

- 0.12mm, 1 x 2 teeth, tip length 2.5mm
- 6.0mm tying platforms
- Colibri style shafts
- Dolphin handle, length 94.0mm

Allows for marking intended incision sites as well as fixation of globe for LRI procedure.

Titanium Instruments

Duckworth & Kent Ltd.

Terence House
7 Marquis Business Centre
Royston Road, Baldock
Herts SG7 6XL England

Tel: +44 (0)1462 893254 Fax: +44 (0)1462 896288
Web: www.duckworth-and-kent.com
Email: info@duckworth-and-kent.com

www.duckworth-and-kent.com

